


# LOCATOR® ATTACHMENTS

*Designed with ease of use  
and incredible durability  
for your patient*



## SELF-ALIGNING FEATURE

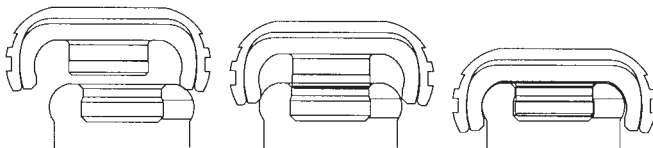


*Self-locating design allows a patient  
to easily seat their overdenture*

### For Root Supported Overdentures

The LOCATOR Root Attachment is a supra-radicular design that has a choice of a straight post, and two angles (10 and 20 degrees) to accommodate divergent roots. Different color-coded retentive males allow for your choice of retention according to the needs of the patient.

The self-correcting alignment of the LOCATOR attachment reduces wear and eliminates damage to the attachment components by directing the proper insertion of the appliance every time (even if the patient bites their denture into place!).

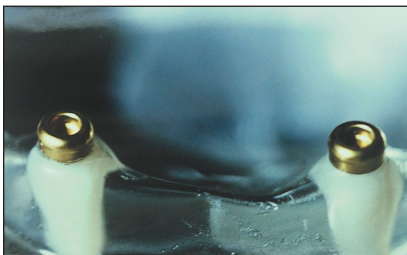

LOCATOR Root Attachments in a model.

### For Use With Cast Copings

A special LOCATOR Cast-To Female is used for cases where a gold coping is desired to cover and protect the root surface. The dental laboratory simply incorporates the stainless steel attachment into a wax-up and casts the gold alloy directly to the LOCATOR attachment.


The LOCATOR has a complete attachment vertical height requirement of only 2.5mm that permits its use in tight inter-arch spaces. In addition, the self-aligning ability of the attachment aids the patient in positioning their prosthesis in a similar manner as a guide plane created by a milled bar. Clinical use has proven the ability of the LOCATOR Attachment to correct patient problems with inserting their overdenture.

To learn more about the LOCATOR  
visit our Web Site at [www.zestanchors.com](http://www.zestanchors.com)

**ZEST ANCHORS, INC.**  
2061 Wineridge Place #100  
Escondido, CA 92029  
Phone: 760-743-7744  
Fax: 760-743-7975


Place orders online or call toll-free  
**1-800-262-2310**

## LONGEST LASTING

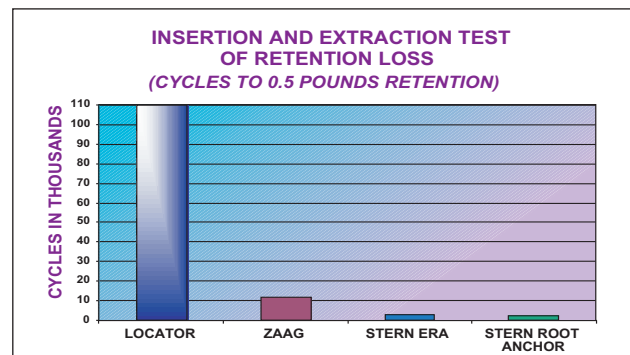

### For Implant Retained Overdentures

The LOCATOR Implant Attachment is designed to reduce the height of abutment restoration on all brands of endosseous implants.

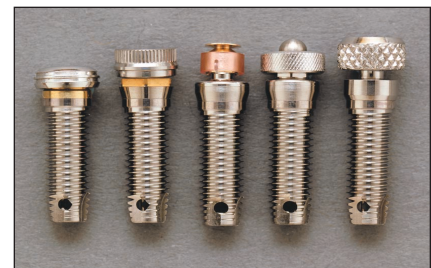

3.17mm LOCATOR    4.85mm ERA    5.82mm Dal-Ro    6.14mm O-Ring    6.22mm EDS

With a total attachment height of only 3.17mm for an externally hexed implant, the LOCATOR saves a minimum of 1.68 mm to 3.05 mm of inter-occlusal space compared to other implant overdenture attachments.

The LOCATOR Implant Attachment can also be used with highly divergent implants. Restoration of a non-parallel implant with angulation of up to 20 degrees (40 degrees of total divergence between two implants) can be accommodated with the LOCATOR system.