


New Product Release – LOCATOR


Introducing the new
**8548 LOCATOR Extra Light
Retention Extended Range Male
(Red) (4 Pack) \$20.00**

New color = RED

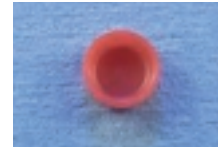

New Item: 8548

Now extra light retention is possible for free-standing implant cases when divergent implants are placed into function and when less retention for multiple attachments or gradual loading of the implants is desired.

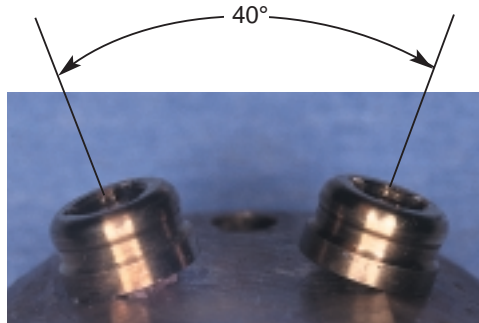

The **new 8548 LOCATOR (extra light retention) Extended Range Replacement Male (Red)** provides less retention (1.5 pounds / 680 grams) and will easily accommodate an off-angled implant (10 to 20 degrees, or 40 degrees between 2 implants).

Our existing 8547 LOCATOR Extended Range Replacement Male (Green) offers more retention for divergent implants (4-3 pounds / 1814-1361 grams) depending upon the degree of angulation. It has less retention with a greater angle.


Like the LOCATOR Extended Range Replacement Male (Green) the **new Red Male** will simply snap into the LOCATOR Denture Cap (metal housing). Plus it has the resiliency and durability that you always expect due to the LOCATOR's stress-free pivoting cap design.

Note: LOCATOR metal housing is available in the LOCATOR Male Processing Package (*Denture Cap with Black Processing Male, Clear & Pink Replacement Males, and Block-Out Spacer*) 8519-2 (2 Pack) or 8519-10 (10 Pack).



Male Processing Pack
8519-2, 8519-10

Extended Range Replacement Males (Red & Green) utilize outside retention only so for implants placed 0-10 degrees please use dual retention nylon males (clear, pink, or blue). LOCATOR Replacement Males are each available separately in a 4 pack for \$20.00.

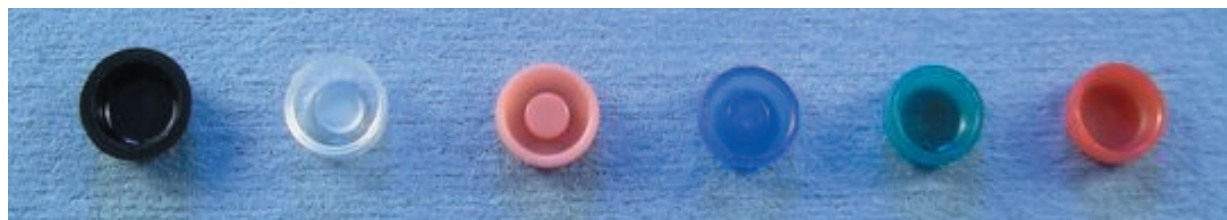

Process
8515

5 lbs.
8524

3 lbs.
8527

1.5 lbs.
8529

4-3 lbs.
8547
(angled)

1.5 lbs.
8548
(angled)

Zest Anchors, Inc. 1.800.262.2310 or 760.743.7744

Place orders online at www.zestanchors.com
and receive free UPS ground shipping.